

# PARA PULSER®

**TAKE BACK YOUR POWER**



## Product description

Ver 4.0

October 2025



# PARAPULSER® – the magnetic-(r)evolution

The **Parapulser®** has been on the market since 2012 and is the most compact, advanced and powerful magnetic pulser or **PEMF** device in its class (**P**ulsed **E**lectro **M**agnetic **F**ield). It is not only used for applications according to Dr. Beck, but also serves to increase physical and mental well-being. Its applications are virtually unlimited: strong magnetic pulses are suitable for activating cell metabolism in all areas of the body. With the **Parapulser®**, the pulse frequency can be programmed to individual frequencies. In addition, the magnetic pulses are charged with information from the quantum field to support the overall effect.

The effect of magnetic pulses is based on the physical effect of magnetic induction, in which all ions (salts dissolved in the body) in the body's cells are set in motion, gently activating the body's own self-healing powers.

The **Parapulser®** is not a medical device and does not replace a visit to a doctor or alternative practitioner.



## Technical highlights:

- **Pulse frequency 1 – 20 Hz:** The **Parapulser®** can generate 1 – 20 magnetic pulses per second. This provides the user an unrivalled range of intensities and effects.
- **Adjustable magnetic pole:** With the **Parapulser®**, the polarity of the magnetic pulses generated can be set to "north", "south" or "alternating pole" at the touch of a button.
- **Individual pulse frequencies:** Up to 12 individual pulse frequencies can be programmed to three digits after the dot in the range 1 – 20 Hz. This gives the user the freedom to program and arrange the most often used pulse frequencies according to their preferences. This reduces the number of button presses required to switch between them.
- **Acoustic pulse counter:** The acoustic pulse counter enables precise and equal dosing of the magnetic pulses over the desired body locations. The acoustic pulse counter can be programmed between 10 and 9999 pulses. Whereby after each counting interval a beep tone is generated.
- **0.66 Tesla pulse strength:** Each magnetic pulse has a strength of 0.66 Tesla (6,600 Gauss) at the output-surface regardless of the pulse frequency.

- **5000 magnetic pulses:** Using state-of-the-art technology, the **Parapulser®** generates up to 5,000 magnetic pulses at a time, even at the highest pulse frequency, before automatically switching off due to self-heating and thus protecting against skin burns.
- **Single button operation:** The **Parapulser®** can be operated with a single button in daily use.

## Display functions:

- **Magnetic pole:** During operation, the polarity (north/south pole) of the generated magnetic pulse is indicated by LEDs (red/blue).
- **Intensity-step:** The currently set intensity-step (1 – 12) is displayed by color LEDs in a coded way.

## Device versions:

The **Parapulser®** is available in three versions:

- 1.) **Parapulser® BASIC (2 Hz)** - for beginners
- 2.) **Parapulser® SEMIPRO (10 Hz)** - for common use
- 3.) **Parapulser® PRO (20 Hz)** - for professional use


The **BASIC** and **SEMIPRO** model offers a *free trial* (\*) to test the enormous power of the **PRO (20 Hz)** version, if **2 Hz** or **10 Hz** are not sufficient to make a change. In addition, the **BASIC** and **SEMIPRO** model can any time be upgraded to a higher version by entering a purchased upgrade-code into the device.

(\*) The free trial includes **40** applications, each with **5,000** magnetic pulses at pulse frequencies up to **20 Hz**.



## Scope of delivery:

The **Parapulser®** comes in an elegant case containing all necessary parts:

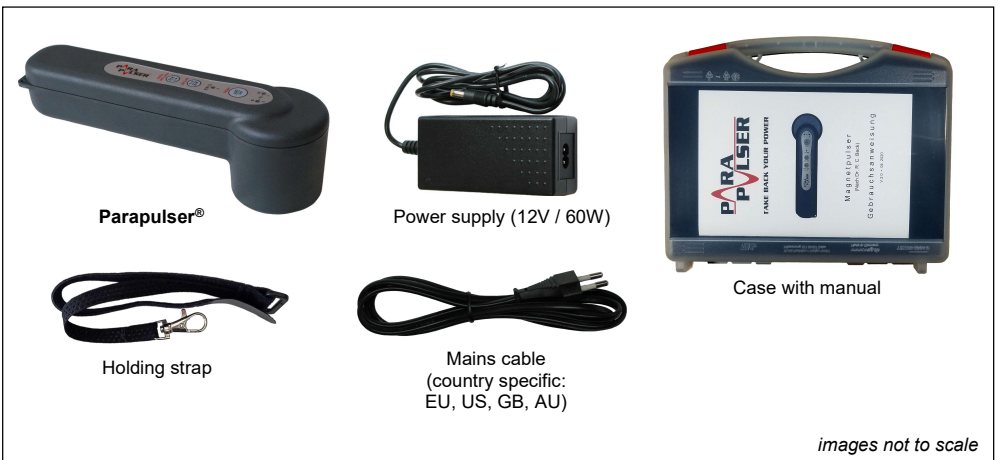

## **Application:**

Boost your well-being with the help of powerful magnetic pulses and regain your energy! The **Parapulser®** supports you more effectively than any other comparable device. It is very easy to use:

- 1.) Switch on the **Parapulser®** by shortly pressing the "on/OFF" button,
- 2.) Start the magnetic pulses by shortly pressing the "on/OFF" button,
- 3.) After use, switch off the **Parapulser®** by pressing the "on/OFF" button for 2 seconds.



Example for application of the **Parapulser®** at the knee joint

---

## **Legal notice:**

In accordance with the German Advertising of Medicines Act (HWG) and the German Unfair Competition Act (UWG), we are obliged to point out that the effectiveness of magnetic field devices is still the subject of controversy in the scientific community and is not yet generally recognised.

Handed out by:

### **Mikas Elektronik Entwicklungen**

Dipl.-Ing. Alexander Mikas

Bahnhofstr. 15

82340 Feldafing

Germany

Phone : +49 8157- 92 42 92

Mail : [info@mikas-elektronik.de](mailto:info@mikas-elektronik.de)

Web : [www.mikas-elektronik.com](http://www.mikas-elektronik.com)

WEEE-Reg.-Nr. DE 13571115